


Cobra VIS OCT Spectrometer Series

Unveil a new dimension in 3D imaging

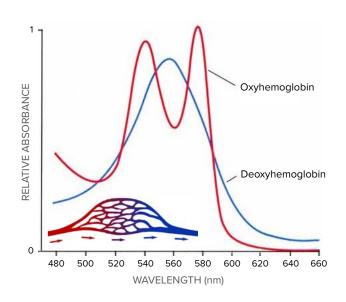

The differential in absorption at 550 nm between oxyhemoglobin and deoxyhemoglobin provides additional information about vascular perfusion.

MAXIMIZE INFORMATION

Collect blood perfusion data and ultra-high resolution OCT images simultaneously at 550 nm

Industry-leading line scan rate: up to 250 kHz

Low crosstalk camera & diffraction-limited optics for exceptional roll-off; as little as 6 dB at 1.2 mm

Proprietary VPH gratings; highly efficient optical design with superior SNR & subpixel resolution

USB 3.0 interface option for 20-130 kHz models allows easy, cost effective system integration

OEM ready: robust, compact & athermal

Access all the benefits of vis-OCT with our newest addition to the Cobra spectrometer family.

The Cobra VIS has been designed with the same industry-leading performance as our flagship Cobra-S spectrometer for high resolution, high speed imaging at visible wavelengths. It gives you access to valuable oxygenation data while nearly doubling resolution as compared to 800 nm OCT imaging, empowering the structural and functional imaging of tissues. We conquer motion artifacts with 250 kHz imaging and offer convenience with a USB 3.0 interface. From 3D blood oxygenation mapping and microscopy to high resolution retinal imaging, the Cobra VIS gives you more information than SS-OCT ever could.

ANGIOGRAPHY | 3D MICROSCOPY | HIGH-RES RETINAL IMAGING




We put high resolution vis-OCT easily within reach.

Map perfusion and tissue oxygenation confidently by using the best possible optical design and camera on the market – diffraction limited and balanced for chromatic aberration. Try our new USB 3.0 interface and dedicated software libraries for high resolution, high-speed imaging from the convenience of a laptop!

	CS550-555/80	CS550-573/115	CS550-600/200	CS550-475/150
Imaging Depth	1.9 mm	1.3 mm	0.9 mm	0.7 mm
Wavelength Range	515-595 nm	515-630 nm	500-700 nm	400-550 nm
Bandwidth	80 nm	115 nm	200 nm	150 nm
Spectral Resolution	0.04 nm	0.06 nm	0.1 nm	0.08 nm
Max Line Rate	20 kHz, 80 kHz, 130 kHz, or 250 kHz			
Pixels	2048			
Interface	USB 3.0 (up to 130 kHz) or Camera Link (up to 250 kHz)			
Dimensions	22 x 12 x 6.5 cm			
Weight	1.3 kg			
				■ LE COSCOSCOSCOSCOSCOS

INCREASING AXIAL RESOLUTION >

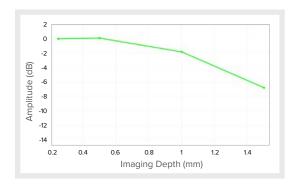

The Cobra VIS uses proprietary optics and a high sensitivity camera to deliver better roll-off than most OCT spectrometers, improving your resolution at increasing depths.

OEM CUSTOMIZATION

At Wasatch Photonics, we go beyond custom to create unique, bespoke spectrometers for our OEMs. We'll share our deep understanding of spectrometer and OCT system design, working with you as a collaborator to create products to differentiate you in your marketplace. From custom cameras, gratings, and athermal lens sets to drop-in replacements for legacy designs, we can develop a solution optimized for your imaging needs.

Contact us for greater clarity

WP-PS_Cobra-VIS_01Apr19

